Installation and operation guide
 Instructions pour l'installation et l'emploi
 Bedienungsanleitung
 Istruzioni per l'installazione e l'uso
 Instrucciones para su instalación y uso
 Instruções para a instalação e o uso
 Installations och användarmanual

Digital sterilization data recorder
 Enregistreur numérique des données de stérilisation
 Digitaler Zyklusspeicher
 Registratore digitale dei dati di sterilizzazione
 Unidad de registro para ciclos de esterilización
 Registrador digital dos ciclos de esterilização
 Digital lagring av cykeldata

Manufacturer:

W&H Sterilization S.r.l
 Italy, I-24066 Pedrengo (BG)
 Via delle Crocette 2
 t +39/035/65 94 111
 f +39/035/66 83 87
 e office@wnhster.it
 wh.com

Importer:

Lisa log⁶⁴

DIGITAL STERILIZATION DATA RECORDER
Installation and operation guide

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Installation and operation guide

Dear Owner,

Thank you for purchasing our digital data recorder **Lisa log 64**.

W&H Sterilization is pleased to provide information concerning the use of this equipment herewith.

We remind you that in order to use the digital recorder correctly, it is necessary to read this manual carefully. All our products comply with the current safety standards and do not involve any risk for the operator when used according to the instructions provided.

Whilst wishing you every success in your work, please note that the reproduction of this manual is forbidden and that, owing to continuous research and technological development, the equipment specifications could be changed without prior notice.

W&H Sterilization S.r.l.

TABLE OF CONTENTS

INTRODUCTION	2	DESCRIPTION AND USE OF THE SYSTEM	7
General Warnings.....	2	LISA log 64 operating principle	7
About the Warranty	2	Use	7
PACKAGE CONTENTS.....	3	Memory card	7
General notes on delivery.....	3	At the start of the sterilization cycle	7
Packing size and weight	3	At the end of the sterilization cycle.....	7
List of the components	3	Data transfer	7
TECHNICAL DATA	4	Signal description	8
Characteristics	4	The files.....	10
LISA sterilizer compatibility	4	Control code.....	10
INSTALLATION	5	File name.....	10
Connections	5	Directory structure	11
Positioning the Modules	6	MAINTENANCE.....	12
LISA sterilizer configuration.....	6	External cleaning	12
Connection of USB reader to PC/MAC computer	6	Spare parts	12

INTRODUCTION

This manual provides instructions for:

- safe and efficient operation of the equipment;
- continuous and scheduled maintenance.

The sizes mentioned in the manual are not binding. Drawings and any other documents delivered with the equipment belong to W&H Sterilization, who reserves all the rights, and they cannot be made available to third parties.

The reproduction, either total or in part, of the text and the illustrations, is forbidden.

General Warnings

- The equipment should be used in compliance with the procedures described in the manual and never for purposes other than those detailed in the manual.
- The user is responsible for installing, operating and servicing the equipment in accordance with the instructions listed in this manual.
- The manufacturer is not liable for any direct and/or indirect damages due to the loss of data.

Indication of use: This manual is an integral part of the product and must be kept close to the equipment for reference. The equipment was designed for professional use. The equipment is to be used only for the purpose for which it was designed:

Lisa log 64: Device designed for the digital recording on memory card (MMC) of test - and sterilization cycles run on LISA autoclaves. The data can be transferred to computers using the appropriate reader (USB-reader included)

About the Warranty

W&H Sterilization guarantees the Lisa log 64 recorder, if used in accordance with the instructions described in this manual, for a period of **12 months** for both mechanical and electrical components. For those components selected and purchased by us (i.e. "MMC" and USB reader) please consider the warranty included in the single packaging.

The warranty starts from the date of delivery to the customer, which must be proven by returning the warranty form duly filled in and signed.

In case of dispute, the date on the delivery note relating to the equipment serial number will be valid. The warranty does not cover defects or damage caused by:

- Incorrect installation
- Inadequate maintenance
- Tampering or unauthorized modifications
- Incorrect use of the equipment
- Non-compliance with the instructions included in this manual

The replacement of equipment or the extension of the warranty period because of breakdown is excluded

Repair under warranty is carried out at W&H Sterilization's discretion and does not include packing and shipment costs.

	GB	WARRANTY	FR	BON DE GARANTIE
	DE	GARANTIE	IT	GARANZIA

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PACKAGE CONTENTS

General notes on delivery

Upon delivery of the equipment check the condition of the package and keep it in case of future transit.

Open the package and check that:


- the contents correspond to the delivery note
- there is no evident damage.

I In the event of damage or missing items, inform the carrier and the dental dealer or W&H Sterilization immediately.

Packing size and weight

Size (L x D x H)	412 x 317 x 70 mm
Weight	1,4Kg

List of the components

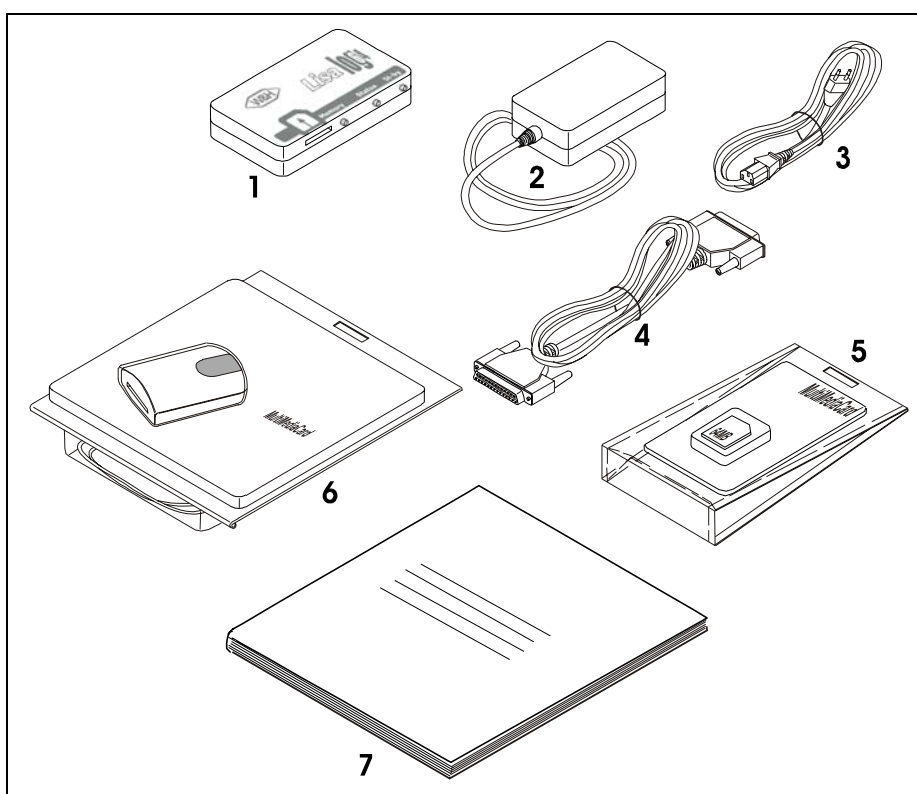Item	Description
1.	Lisa log 64 recorder
2.	Power supply
3.	Mains cable
4.	Parallel interface cable
5.	Memory MultiMediaCard Kit (*)
6.	MultiMediaCard USB reader Kit (*)
7.	Installation and operation guide (this manual)

(*) The memory card and the USB reader of the supply are properly selected by W&H for the use with LisaLog64 recorder. Please refer to the instructions (hard copy or CD version from the corresponding supplier (s)).

TECHNICAL DATA

Characteristics

Lisa log 64 recorder

Operating temperature	+0°C ÷ +40°C (avoiding heat sources and humidity)
Storage temperature	-20°C ÷ +40°C (avoiding heat sources and humidity)
Sizes (without memory card and interface cables)	126 x 68,7 x 24,6 mm (LxDxH)
Weight (without memory card/cables)	120g
Use	internal
Power supply voltage	12Vcc (by the supplied power supply module)
Max power consumption	1,2W

POWER SUPPLY

Operating temperature	+0°C ÷ +40°C
Input voltage	110~240Vac 50/60Hz
Output voltage	12Vcc +/-10%
Output current	Min. 200mA

MEMORY CARD

Technology	MMC
Capacity	64Mb
Formatting	FAT16 – FAT32 preformatted

MINIMUM PC/MAC HARDWARE REQUIREMENTS

PC/MAC	1 available port USB1.1 (or higher)
PC/MAC	CD-ROM reader
PC operating system	WINDOWS 98SE, 2000, ME, XP or later
MAC operating system	OS 9.1.x, OS X v10.1.2+

Brand and products are trademarks of their respective holder/s

LISA sterilizer compatibility

Lisa log 64 is designed for the use with all LISA MB17 and LISA MB22 sterilizers on the market, and does not require any software upgrade of the sterilizer.

INSTALLATION


Caution: Before starting the installation, turn off the sterilizer and disconnect the mains cable from the socket.

Connections

Install **Lisa log 64** digital recorder as follows:

- Turn off the LISA sterilizer;
- Connect the parallel interface cable between the printer port of the sterilizer and the parallel port of the recorder module;
- Plug the power supply output cable in the **Lisa log 64** socket;
- Connect the mains cable between the power supply module and an available mains socket;
- Switch on the sterilizer;
- Plug the memory MultiMediaCard into the slot of the **Lisa log 64** in the orientation indicated in the diagram.

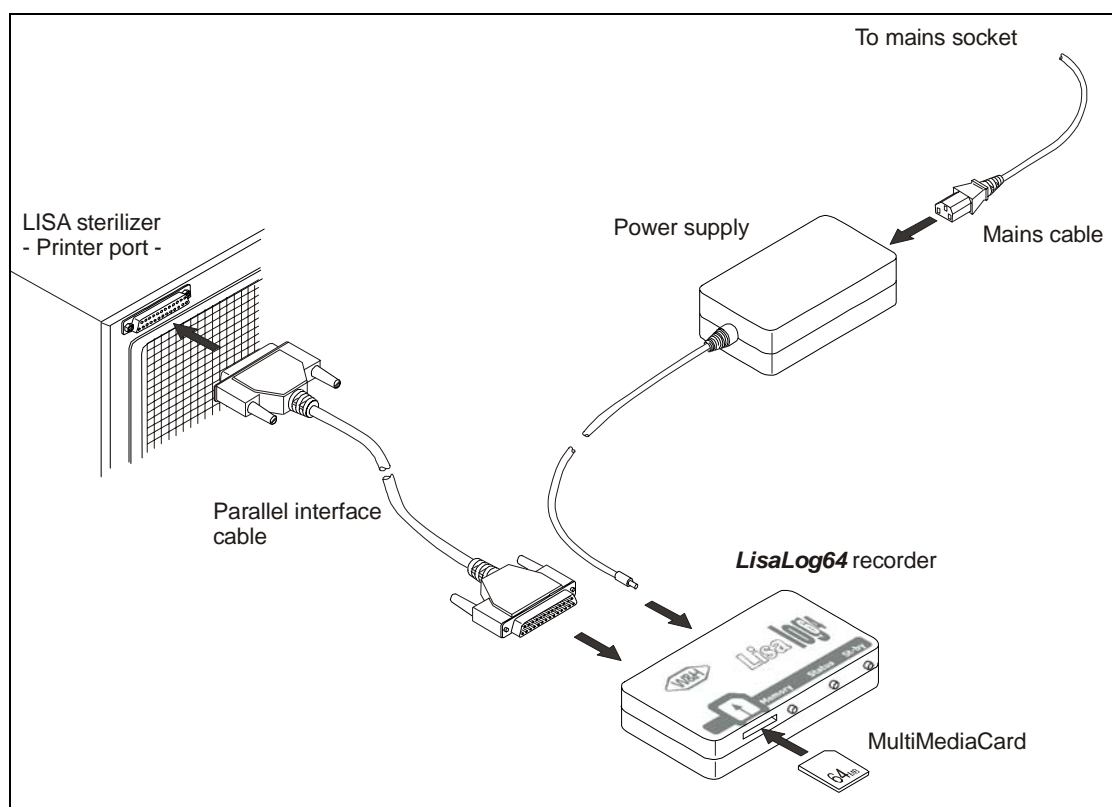Positioning the Modules


CAUTION:

Do not arrange the modules near heat sources (i.e. behind the sterilizer).
On the sterilizer do not obstruct the aeration holes of the water recovery tank.


CAUTION:

Install **Lisa log 64** so that the signalling is easily visible and abnormal conditions promptly identified and possible data loss avoided (e.g. power supply black out or incomplete data saving).

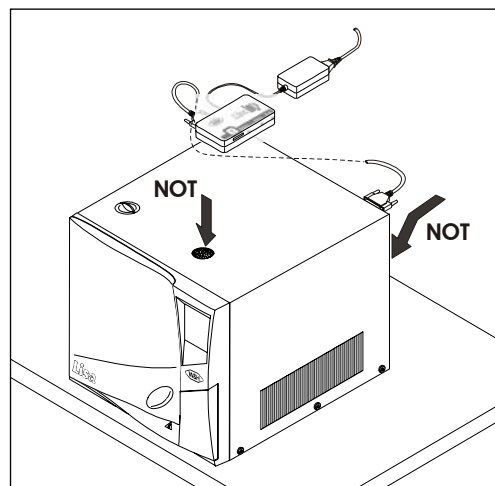Lisa log 64 can be installed and can operate correctly in any position; we suggest however not to use the module with the memory card towards the bottom in order to avoid unplugging the card through vibrations or impacts.

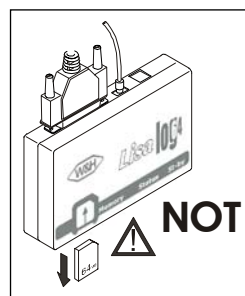LISA sterilizer configuration

No setting required.

Lisa log 64 is compatible with all LISA MB17 and LISA MB22 sterilizer. Moreover which printer is selected in the sterilizer software is not relevant.

Connection of USB reader to PC/MAC computer

The saved data on the MMC memory can be read, copied and printed with a PC or Macintosh computer equipped with a compatible interface slot.

The USB reader supplied allows the use of MMC memories on PC/MAC not provided with MMC reading slot and equipped with a port USB 1,1 or higher.

To install the USB reader on your PC/MAC, proceed as follows:

- Turn on the PC/MAC;
- Connect the USB cable between the reader and the available USB port of the PC/MAC;
- Follow the instructions of the operating system.

N.B. *Most of the operating systems recognize the reader and automatically proceed with the installation of the device. Other operating systems request the installation of the disc; in these cases follow the instructions to complete the installation; for any additional information refer to the handbook supplied with the USB reader.*

Upon completion of the installation of the USB reader driver, the MMC card will be seen as a removable peripheral memory device on which it is possible to perform any file copy, movement, reorganization and opening operations.

CAUTION:

Before unplugging the memory from the USB reader, make sure that all the operations being performed have been ended, otherwise it is possible to damage the memory or lose the data. Some operating systems allow for the removal of the peripheral device safely via a proper function.

DESCRIPTION AND USE OF THE SYSTEM

LISA log 64 operating principle

Lisa log 64 is a device designed for the use with LISA MB17 and LISA MB22 sterilizers. When connected to the parallel printer port, it saves one file for each sterilization cycle run.

Use

In normal conditions the **Lisa log 64** does not require any intervention by the user; should the device detect any abnormal conditions or possible loss of data (cable detached, memory card full or failed, etc.), the recorder provides proper visual and/or acoustic signalling as detailed in the next section.

At the end, the device switches automatically into sleep-mode after 5 minutes of inactivity. The exit from stand-by will occur automatically as soon as new data is recorded, the memory card is removed or any activity is detected.

Memory card

The memory card provided is already formatted according to FAT16 or FAT32 standard.

When the memory card is plugged in, the **Lisa log 64** performs a check before signalling that the card is valid.

WARNING:

The time necessary for the memory check ranges from 2 s and 40 s or more depending on the memory type and formatting type.

CAUTION:

Do not remove or plug in the memory card during a sterilization cycle.

Lisa log 64 is able to memorize a complete sterilization cycle without the memory card plugged in. To save the sterilization cycle on the memory card the card must be inserted after the sterilization cycle.

At the start of the sterilization cycle

At the beginning of the sterilization cycle, (once preheated) the Memory LED flashes green. It means that **Lisa log 64** is receiving data from the sterilizer.

At the end of the sterilization cycle


At the end of the sterilization cycle, the Memory LED remains green. The recorder switches into sleep-mode after 5 minutes.

WARNING:

If at the end of the sterilization cycle the Status LED flashes orange, it means that some data was incomplete or missing.

In order not to lose the previous cycle, check the connection of the parallel cable, and prior to opening the sterilizer door perform the re-printing as described below:

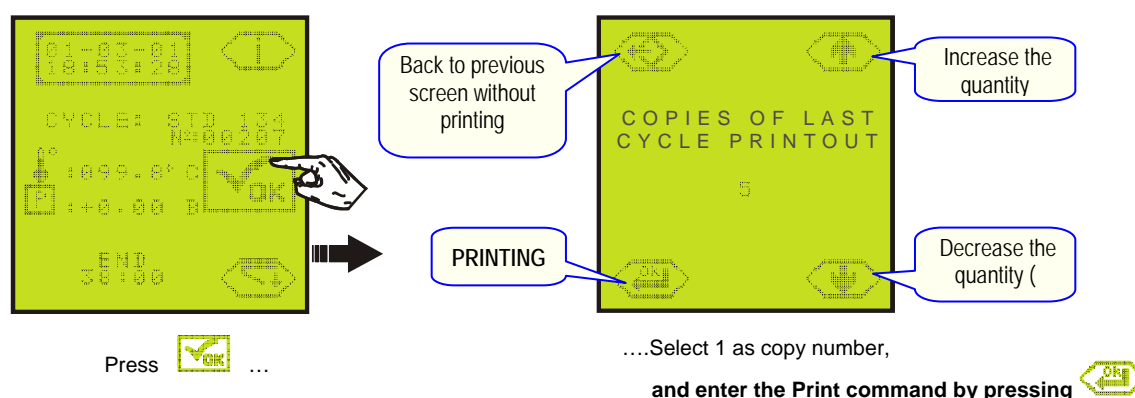The STATUS LED will stop flashing orange when a new complete cycle is saved or the last one is reprinted correctly as indicated above. If the condition remains even after the reprinting function, refer to the Signaling table.

Data transfer

We recommend that the user transfers cycle data from the card to the computer on a regular basis.


Signal description

LED

Three Led's are available for a quick view of the recorder condition:

- **Memory** – Status LED of the memory card,
- **Status** - Status LED of the recorder module,
- **St-by** - Stand-by status LED

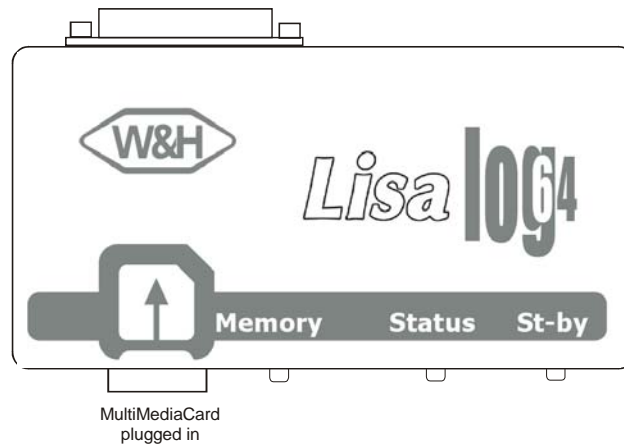The presence of acoustic and/or red signalling means that there are abnormalities and/or possible loss of data. Other combinations of signalling are not associated with anomalies, but provide an overview of the operating status of the device.

Led			MEANING
Memory	Status	St-by	
Off	Off	Off	Memory card missing, Stand-by status or power supply missing
Red blinking + MMC missing			Memory card missing during the cycle - To save the data, plug in the card as soon as the sterilization cycle ends. WARNING:  Do not remove or plug in the memory card during a sterilization cycle. <i>Lisa log 64</i> is able to memorize a complete sterilization cycle without memory card plugged in. To save the sterilization cycle on the memory card the card must be inserted after the sterilization cycle.
Green steady			It is possible to start a sterilization cycle, and/or transfer the data onto the computer
Green blinking			<i>Lisa log 64</i> is receiving data. Do not unplug the memory card. The data will be saved onto the memory card at the end of the cycle, indicated with a faster blinking of the LED. Note: If the cycle is not ended and no data is being received for more than 45 minutes (Sterilizer Off) the file will be saved as incomplete.
Red/Green blinking			The diagnostic of the memory card is in process (from 2 to 40 seconds)
Red blinking + MMC plugged in			The memory card is defective or not formatted
Red steady			The memory card is full. It is impossible to save other cycle data! WARNING:  Do not remove or plug in the memory card during a sterilization cycle. Even if the memory card is full, <i>Lisa log 64</i> is able to memorize a complete sterilization cycle. Insert an empty card only at the end of the sterilization cycle.

Led			MEANING
Memory	Status	St-by	
Off	Off	Off	Stand-by status or power supply is missing.
	Green steady		OK – the data can be saved on the memory card
	Red/Green blinking		OK – the data can be saved, but the remaining memory capacity is less than 50 cycles. Transfer the data onto your PC and insert a new one. It is also possible to erase data from the previous one and re-use it.
	Red blinking		Impossible to save data – cable disconnected.
	Red steady		Impossible to save data– see Memory LED status.
	Orange blinking		Warning: the last file saved does not contain a complete data cycle. Check the parallel cable and perform a reprint (see page 7). Note: <i>The signalling continues until a complete file is saved; during a cycle this signalling does not mean that there is an abnormality.</i>
	Off	Off	With at least one LED on, the recorder is OK. With all the LEDs off, no power supply
		Yellow	OK – Stand-by status
Orange	Orange	Off	Power supply problems. Check the cable connections (Connect and disconnect). If the problem remains, Lisa log 64 or the power supply may be damaged – call for service. Note: <i>When Lisa log 64 and a printer (via a special adaptor) are used at the same time and the printer is turned on, the signalling may remain on even if the power supply is not connected to the Lisa log 64.</i>
Orange	Orange	Yellow	Device not functioning (try to turn it off and then back on)

ACOUSTIC SIGNALING

Signaling	Description
4 beeps every 15 minutes	CABLE DISCONNECTED The parallel cable is disconnected. If a cycle is launched the data will not be saved. Check the connection between the sterilizer and Lisa log 64 .
3 beeps at the beginning of the cycle	MEMORY NEARLY FULL It is possible to save data, but the remaining memory capacity is less than 50 cycles. Transfer the data onto your PC and insert a new one. It is also possible to erase data from the previous one and re-use it.
3 beeps each time the data is received	MEMORY FULL, MISSING or FAILED Impossible to save data. Even if the memory card is full, Lisa log 64 is able to memorize a complete sterilization cycle. Insert an empty card only at the end of the sterilization cycle. <div style="display: flex; align-items: center;"> <div> CAUTION: If other cycles are started the new received data could be lost. Transfer the data onto your PC and insert a new one. It is also possible to erase data from the previous one and re-use it. </div> </div>

The files

The sterilization data is saved on files with .TXT format and can be opened, viewed and printed by means of a standard text editor. The files are saved with a read-only attribute.

Control code

At the end of the file a control code is included in order to allow the user to check the authenticity of the content. The control code can be checked by W&H Sterilization. The control code contains some additional information.

Example:

W&H O.D.Code 00-EC6FF-F62333696D4AD51A

W&H O.D.Code	0	0	-	EC6FF	-	F62333696D4AD51A
1	2	3		4		5

- 1: W&H Original Data Code
- 2: **Lisa log 64** internal programme
- 3: Upgrade loader release
- 4: Resuming code of the **Lisa log 64** status (for W&H service)
- 5: Control code for the file content authenticity check

File name

The file names summarize the main information contained in the file according to the following structure. File name example:

040128 28-01-04 08_41_22 A99 22E02-0095 01898 B&D _ HELIX TEST .txt

The file name is made up of the following information:

040128	28-01-04	08_41_22	A99	22	E	02-0095	01898	B&D _ HELIX TEST	.txt
1	2	3	4	5	6	7	8	9	10

- 1: Date of the cycle: format YYMMDD in order to save the files chronologically:

YY last digits of the year
MM month
DD day

- 2: Date of the cycle: same format used by the sterilizer:

DD-MM-YY for EU versions
YY-MM-DD for JP versions
MM-DD-YY for US versions

- 3: Time of the cycle: same format used by the sterilizer:

HH-MM-SS for EU and JP versions
HH-MM-AM/PM for US versions

- 4: Alarm: (only in case of report with alarm)

with format AXX, where:
XX identifies the alarm code detected by the sterilizer

- 5: Model of the sterilizer:

17 LISA MB17 sterilizer
22 LISA MB22 sterilizer

- 6: Sterilizer programme release

E EU version
U US version
J JP version

- 7: Sterilizer Serial Number

- 8: Sterilizer cycle number

- 9: Cycle type

- 10: File name extension

W&H STERILIZATION				5	6
LISA MB				22	EU 7.1
Dr.				7	
SN. :	02-0095				

Data :	28-01-04 08:13:57				
Cycle :	B&D / HELIX TEST 9				
Number :	01898 8				

PHASE	TIME	Temp C	P	Bar	
Start	00:00	050.0		-0.03	
PV1	05:25 05:25	050.0		-0.96	
PP1	05:28 00:03	050.0		-0.67	
.....					
LEV	27:24 00:30	104.1		-0.03	
.....					
END	27:24	104.1		0.03	
1 2 3					

Data :	28-01-04 08:41:22				
STERIL LOAD					
=====					
W&H O.D.Code 00-EC6FF-F62333696D4AD51A					

PV2	06:39 01:11	057.6		-0.85	
PP2	06:59 00:20	079.4		+0.07	
PV3	08:01 01:02	069.3		-0.75	
PP3	08:05 00:04	080.1		-0.43	
4					

A99	PHASE: PP3				

WARNING:

The file name does not include forbidden characters \ / : * ? " < > | as they violate the rules for the file name generation. All forbidden characters are replaced by the character "_" (underscore). For example, the file name **B&D / HELIX** is converted to **B&D _ HELIX**.

All data forming the file name are obtained directly from the print report and maintain the language used on the sterilizer.

Directory structure

TYPICAL STRUCTURE:


The files are saved on the MMC card according to the following directory structure:


```

Lisa + serial number
├── Year (four digits)
│   └── Month (two digits)

```

example of a typical structure automatically created for a sterilizer:

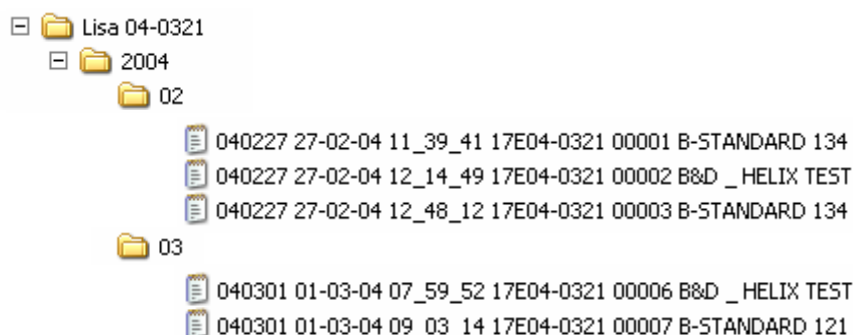
ADDITIONAL FOLDERS IN CASE OF ALARMS:

In case of alarms the folder Axx will be created inside the general folder of the sterilizer generating the alarm; the file will be saved here with the alarm code:

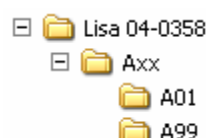The incomplete files (from which at least the sterilizer's serial number can be obtained) will be placed in a folder named "xxx":

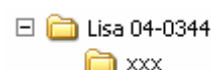ADDITIONAL FOLDERS IN CASE OF INCOMPLETE FILES:

The reception of incomplete data (e.g. due to cable disconnection or power supply blackout) could render the complete construction of the file name according to the structure described impossible. If so the file name will use the data of the previous cycle instead of the lacking information and will show in brackets the incomplete information: such expedients enable the period in which this data has been received to be found.

Example:

If after a file with file name "040128 28-01-04 08_41_22 22E02-0095 01898 B&D _ HELIX TEST .txt" incomplete information is received and where only date and time are identifiable (see example)

DVe	26:14	04:00	107.2	-0.88
VEN	26:54	00:40	106.0	-0.11
LEV	27:24	00:30	104.1	-0.03


END	27:24		104.1	-0.03

Date :	28-01-04 11:24:37			
	STERIL LOAD			
	=====			


W&H O.D.Code 00-EC630-F62333696D4AD51A				

the file name will consist of date and time of the report and in brackets some information of the previous file "040128 28-01-04 11_24_37 (22E02-0095 01898 B&D _ HELIX TEST).txt"

The incomplete files and/or alarms from which the sterilizer's serial number cannot be obtained will be placed in a general folder with serial number "xx-xxxx":

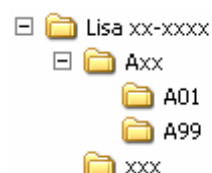An example of the directory structure including various cases previously described:






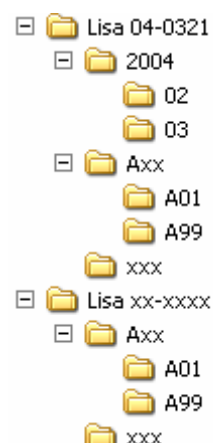
MAINTENANCE

External cleaning

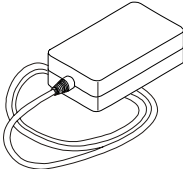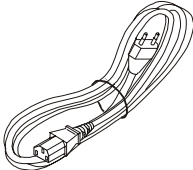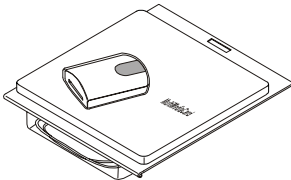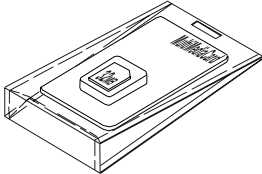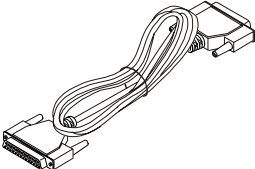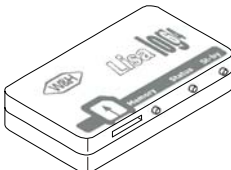Clean with a damp cloth that is not abrasive and does not contain corrosive materials (pH neutral).

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Vänligen notera att denna manual är copyrightskyddad och att på grund av ständigt pågående utvecklingsarbete kan produktens specifikation komma att ändras utan att meddelande ges innan.